GERMAN
POOL

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手部按摩器
Hand Massager

MP-60

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WARNING

Important

For safety purposes, you should read the instructions carefully before using the appliance. Retain the instructions for future reference.

The Hand Massager is designed for domestic use and may only be operated in accordance with these instructions.

1. DO NOT use the massager in extreme heat, high humidity, or near an open flame.
2. DO NOT use the massager if you suffer indisposed symptoms such as inflammation and/or blood stasis. DO NOT use the massager if you have had surgical operations.
3. Remove all hand jewellery before using the massager.
4. DO NOT apply nail polish an hour before using this massager.
5. Place your hand inside the massager before switching it on. This is to avoid damage to the air bag due to excess pressure. It is advisable that users wear disposable plastic gloves whilst using this machine.
6. Please use this massager's supplied adapter and observe basic safety precautions when using electricity.
7. Take out the battery when it has run out. This is to avoid damage caused by battery fluid leakage. Dispose of used batteries responsibly.
8. The hand massager is a household product; it is not suitable for diagnostic, therapeutic, or commercial purposes.
9. DO NOT use aqueous cleaners to clean the device. To clean, wipe the product with a sterile cotton wool.
10. The product has heat massage modes. Users insensitive to temperature differences should take care when using this product.
11. DO NOT use this product while driving.
12. Remove the battery when the massager idles for long periods without use to avoid damage caused by battery fluid leakage.

INTRODUCTION

The MP-60 Hand Massager's intelligent air pressure technology is controlled by a built-in microchip. Through kneading of palm points, massaging motions, and Far Infrared Rays (FIR) heat compression, the product helps promote blood circulation, accelerate metabolism, and soothe aches or pains in the finger joints, tendons, and muscles.

By kneading special pressure points on the palm, the massager can relieve tension and stress, combat fatigue, promote blood circulation, and help regulate the function of internal organs, keeping an optimum bodily function.

The product is suitable for:

- Office workers who are confined to a desk, and those who operate computers for long periods of time
- Artists who frequently use their fingers (e.g. pianists)
- Anyone who suffers from aching palms or hand fatigue
- Anyone who experiences general fatigue
- Anyone who lives a largely sedentary lifestyle

PACKAGE CONTENT


The contents of this package include the Main Unit, a Power Adapter, batteries, a user manual, a carrying pouch, and a pair of plastic gloves.

Main Unit


User Manual

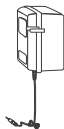Power Adapter

Batteries (AAx4)

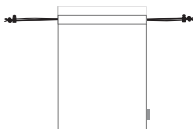Carrying Pouch


Plastic Gloves

* Specifications are subject to change without prior notice. Refer to www.germanpool.com for the most up-to-date version.

PRODUCT STRUCTURE

Main Unit

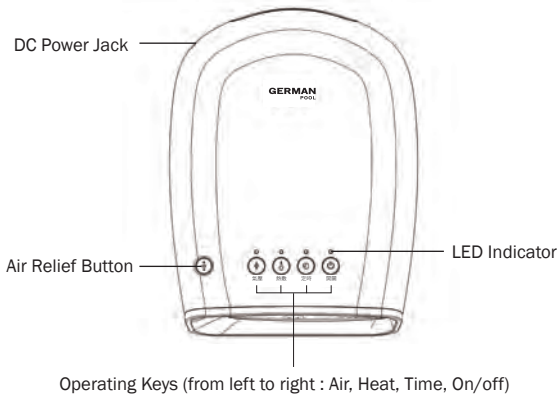* Specifications are subject to change without prior notice. Refer to www.germanpool.com for the most up-to-date version.

OPERATING INSTRUCTIONS

Please check the integrity of the Main Unit, the Power Adapter, and supplied accessories before using the product.

Choose from one of the two following options to start the product:

1. Batteries

- a. Open the battery cover on the back of the Main Unit.
- b. Install four AA batteries. Pay attention to their polarities.
- c. Close the battery cover.

2. Power Adapter

- a. Insert the cylindrical connector of the Power Adapter into the Main Unit's DC Power Jack.
- b. Connect the other end to an AC power socket. The LED will switch on once power is connected.

Starting Up

Press and hold the  key for 5-6 seconds.

The default settings upon start-up are as follows:

- Air: Circulation
- Heat: Off
- Time: 5 minutes

To turn off the Main Unit, press the  key.

To manually adjust the Main Unit, press the related function key to switch air pressure and FIR heat modes.

OPERATING INSTRUCTIONS

Manual Adjustment

Press these function keys to switch between Air Pressure Mode, FIR Heating Mode and timing options:

- ⊕ to select between four pressure settings
- ⊕ to turn on FIR heating and select between two heating options
- ⊕ to select between 5, 10, and 15-minute operation cycles

It is possible to adjust the above at any time during operation. The LED Indicator lights will flash a corresponding number of times to indicate the chosen mode:

Air Pressure Mode ⊕	FIR Heating Mode ⊕	Time ⊕
Circulation (1 flash)	Off (off)	5min (1 flash)
Soft (2 flash)	Low Heat (on)	10min (2 flash)
Standard (3 flash)	High Heat (2 flash)	15min (3 flash)
Powerful (4 flash)		

Note: Press the Air Relief Button ⊕ to decompress if the pressure is too high.

CARE & MAINTENANCE

Storage

Place the unit in a safe place after use. In case the storage period exceeds 6 months, please turn on the Main Unit to check its function. Do not store close to wall, floor, or roof. Keep in good ventilation conditions, away from direct sunlight and/or corrosive substances. Take measures to prevent damage to the unit from moisture, rats, and pollution. Do not put heavy objects on the device.

Transportation

The product has a simple shockproof design and can be transported by air, railway, highway or ship. Avoid rain, snow, heavy impacts, and always keep the unit the right side up.

Check the power supply

Check the power supply before use. Do not switch on the device if the voltage exceeds the rated range of the device.

Maintenance

Contact a qualified technician if there is any fault with the unit. Do not disassemble or attempt to repair the device yourself.

TROUBLESHOOTING

If there is any problem when you use the device, please follow the instructions below to check. If the problem persists or the device does not work normally, please turn off the power supply and contact service centre or a licensed technician. Do not attempt to disassemble and repair the device yourself.

Fault	Possible Cause	Solution
Unit fails to power on	Battery is low	Change batteries
	Batteries are not correctly installed	Ensure batteries are correctly placed in accordance with the positive and negative polarity symbols displayed inside the battery compartment
	Attempting to start up during the power off sequence	The product will start up after it completes the power off sequence
	Power Adapter is not properly connected	Ensure it is properly connected to the main unit
Air pressure or heat compression is low	Battery is low	Change batteries Ensure batteries are correctly placed in the battery compartment

SPECIFICATIONS

Model	MP60
Voltage	6V
Power	6W
Ideal Environment	Temperature : 5°C-40°C
	Relative Humidity : ≤80%
	Atmospheric Pressure : 860hpa~1060hpa
Dimensions	(D) 265mm
	(W) 201mm
	(H) 90mm
Net Weight	1129g

Specifications are subject to change without prior notice. Refer to www.germanpool.com for the most up-to-date version.

If there is any inconsistency ambiguity between the English version and Chinese version, the Chinese version shall prevail.

請填寫以下空格並與發票一併保存，以便維修時核對資料。用戶提供本公司之所有資料必須正確無誤，否則免費保用將會無效。
Please fill out the form below. This information and the purchase invoice will be required for any repairs. Warranty will be invalid if information provided is found to be inaccurate.

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Warranty Terms & Conditions

Under normal operations, the product will receive a 1-year full warranty provided by German Pool (effective starting from the date of purchase).

* This warranty is not valid until Customer's registration information is returned to and received by our Service Centre within 10 days of purchase via one of the following methods:

- 1) Visit our website and register online:
www.germanpool.com/warranty.
- 2) Complete all fields on the attached Warranty Card and mail it back to our Customer Service Centre along with a photocopy of the purchase invoice (Warranty Card and invoice copy will not be returned to Customer).

1. Customer who fails to present original purchase invoice will be disqualified for free warranty.
2. Customer should always follow the operating instructions. This warranty does not apply to:
 - labour costs for on-site installation, check-up, repair, replacement of parts and other transportation costs;
 - damages caused by accidents of any kind (including transport and others);
 - operating failures resulted from incorrect voltage, improper operations and unauthorized installations or repairs.
3. This warranty is invalid if:
 - purchase invoice is modified by unauthorized party;
 - product is used for any commercial or industrial applications;
 - product is repaired or modified by unauthorized personnel, or unauthorized parts are installed;
 - serial number is modified, damaged or removed from the product.
4. German Pool will, in its discretion, repair or replace any defective part.
5. This warranty will be void if there is any transfer of ownership from the original purchaser.

Customer Service & Repair Centre

German Pool (Hong Kong) Limited
8/F, China Travel (Cargo) Logistics Centre
1 Cheong Tung Road
Hung Hom, Kowloon, Hong Kong
Tel : +852 2333 6249
Fax : +852 2356 9798
Email : repairs@germanpool.com

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© All rights reserved. Copying, reproducing or using the contents of this manual is not allowed without prior authorization from German Pool, violators will be prosecuted.

www.germanpool.com