GERMAN
POOL

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mini 迷你手提按摩器
massager | mm³

使用說明書 OPERATING INSTRUCTIONS

目錄 TABLE OF CONTENTS

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Please read these instructions and warranty information carefully before use and keep them handy for future reference.

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This massager is a new generation product developed based on physics, bionics and Chinese traditional medicine and clinical practice. Incorporating quintessences of traditional Chinese medicine and adopting vibration and massage functions, it keeps meridians and collaterals clear, activates Qi, improves blood circulation and promotes longevity. It helps to ease chronic pain, acute pain, muscle ache and palsy, and can promote blood circulation, exudation of pain-restricting substances and immediate interdiction of pain transmission.

Scope of Application

It has remarkable soothing effect on chronic pain, acute pain, muscle ache and palsy.

Target Users

Students, computer users, accountants, designers, clerks, teachers and long-hour working personnel at all ages.

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

This device is designed according to the curves of human body. It is equipped with 3 massage heads and a high performance motor achieving the best massage force. The mini motor controlled by computer IC CMOS chip provides the desirable massage effect.

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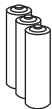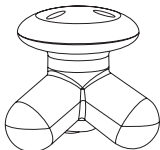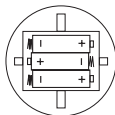
Battery Placement

1. Open main unit cover in anti-clockwise direction. Place three pieces of AAA dry batteries in.
2. Place the battery cover back.
 - Please check if poles of the batteries are correct and if battery cover is securely closed.
 - Please refer to the above figure for positive and negative poles of the batteries.

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Operating Method

1. Place the three kneading massage heads of the massager close to the body part to be kneaded, press the top of the main unit with proper force. The machine starts working.
2. Take the main unit away from the body and give no more force to top of the main unit, the machine stops working.

維修與保養 CARE & MAINTAINACE

- 如本機變髒，可用軟布沾上稀釋酒精擦拭機身。
- 使用後應將按摩器放置於安全地方，切勿在按摩器上放置重物。
- 採取防潮、防鼠、防污染措施。
- 當出現故障時應找專業技術人員修理，切勿隨意拆裝按摩器。
- If the unit body is dirty, clean it with a soft cloth dipped with diluted rubbing alcohol.
- Keep it in a safe place after use. It is strictly prohibited to place heavy articles on top of the device.
- Measures to prevent dampness, mice and pollution are required.
- In case of failure, please have professional technicians repair it. Do not disassemble it by yourself.

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The massager must be used under the following conditions:

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Relative humidity: Max. 80%, without corrosive gas and with good ventilation

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Portable Usage

The product is equipped with simple vibration-proof measures and suitable for transportation by air, railway and ship. Keep it away from rain, snow, splashing and collision. Keep it upright.

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Storage

- If storage period is over 6 months, remove batteries from the unit.
- Do not store it close to ground, wall and roof.
- Good ventilation is required. Keep it away from sunshine and corrosive gas.

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When you run into any problem operating this product, use the table below to troubleshoot the problem. If the problem persists, contact our Repair Centre.

故障現象 Phenomenon	原因分析 Possible Cause	處理方法 Solution
按摩器不啟動。 The massager does not work.	電池是否已消耗盡？ 電池的正負極是否正確？ Batteries using out? Correct poles of the batteries?	更換電池，確認正負極方向，正確插入。 Replace the batteries with correct poles.
按摩力度弱。 Weak massage force.	電池是否消耗盡了？ Batteries using out?	更換電池。 Replace the batteries with correct poles.

- * 以上為常見故障的判斷與檢查，只供用戶作參考之用。
- * 如遇故障切勿自行修理，避免發生危險及對機件造成損害，請立即找專業技術人員或聯絡本公司作出適當維修。
- * Above table enlists the most common problems regarding the operations of the product and their possible causes as well as suggested solutions. It is strictly for users' reference only.
- * Should there be any product malfunction or breakdown, DO NOT repair it yourself. Contact our Repair Centre or consult a licensed technician immediately.

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產品規格及設計如有變更，恕不另行通知。
Specifications are subject to change without prior notice.

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If there is any inconsistency or ambiguity between the English version and the Chinese version, the Chinese version shall prevail.

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Our Repairs & Customer Service Centre will provide rectification of any defect of manufacturing materials free of charge given that such defect is reported within 12 months from date of sale.

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* This warranty is not valid until Customer's registration information is returned to and received by our Service Centre within 10 days of purchase via one of the following methods:

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- 1) Visit our website and register online: www.germanpool.com/warranty
- 2) Mail back a photocopy of the Warranty Registration Card (on the right page) with all fields completed along with a photocopy of the purchase invoice to our Customer Service Centre.

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- Customer should always follow the operating instructions. This warranty does not apply to: damages caused by incorrect voltage, improper operations, unauthorized installation or repairs as well as accidents.
- Customer who fails to present original purchase invoice will be disqualified for free repair services.

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Please mail back a photocopy of this Warranty Registration Card with all fields completed along with a photocopy of the purchase invoice to our Customer Service Centre. Warranty will be invalid if information provided is found to be inaccurate. (Please keep the original copy of this Warranty Registration Card and the purchase invoice, as they are required for any repairs.)

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地址 Address: _____

本人欲定期獲得
I would like to receive:

- 產品資訊
Product News and info
- 折扣優惠
Special Offer

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Repair Centre

German Pool (Hong Kong) Limited
8/F, China Travel (Cargo) Logistics Centre,
1 Cheong Tung Road,
Hung Hom, Kowloon, Hong Kong
Tel: +852 2215 4582
Fax: +852 2356 9798
Email: repairs@germanpool.com

此保用只適用於香港 This warranty is valid only in Hong Kong

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