
GERMAN
POOL

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頸部按摩器
Neck Massager

MN-51

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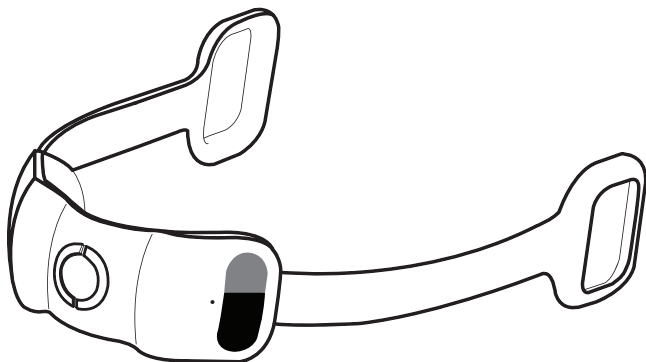READ BEFORE USE

IMPORTANT

For safety purposes, you should read the instructions carefully before using the appliance. Retain the instructions for future reference.

The appliance is designed for domestic use and may only be operated in accordance with these instructions.

1. Do not use the neck massager while driving.
 2. Do not use the massager in a place with high temperature or high humidity.
 3. Do not use this product if you have swelling, inflammation, bruises or any other unusual conditions in your neck. Do not use if you have undergone any neck surgery.
 4. Use this product with care if you suffer from and of the following: severe stroke, hypotension or mental illness. Seek professional medical opinion before use. Stop using immediately if you suffer any discomfort in the process.
 5. Use this product with care and under supervision if you suffer from any physical or mental disabilities, bodily or mental malfunctioning. Stop using immediately if you suffer any discomfort in the process.
 6. Refrain from using the neck massager whilst consuming alcoholic drinks or experiencing sickness.
 7. Please use the power charger provided with this product, and pay attention to the safe use of electricity.
 8. The Neck Massager is designed as a household-use electrical personal care product. It cannot be used for diagnostic, therapeutic or commercial purposes.
 9. Please wipe the massager with sterile cotton wool. Avoid cleaning using aqueous liquids.
 10. This massager emits heat. Take care when using this product if you are insensitive to temperature differences. It is normal for the Neck Massager to generate noise during operation.
 11. The intelligent control on this Neck Massager automatically switches off upon detecting excessive massage pressure. This is normal.
 12. Please tie up long hair before using this Neck Massager.
 13. Please wrap the massager around your neck before switching on.
 14. The massager is a personal care product; it cannot replace the advice and treatment of health care professionals; in case of any questions or need, please contact your health care professionals directly.
- 
- 

INTRODUCTION

The MN-51 Neck Massager uses kneading motions, warm massages, and infrared radiation to relax the pressure points on the neck. Using shiatsu massage techniques and infrared heating, MN-51 helps ease tense muscles, relieve stress, and relax the neck.

The Neck Massager provides fast, effective relief for neck pains, and helps combat symptoms of stress caused by poor blood circulation - such as dizziness or headaches.

The product is suitable for:

- People with neck or shoulder stiffness, or muscle fatigue
- People with neck or shoulder pain, who suffer from numb hands
- People who suffer from headaches or dizziness due to poor blood circulation to the brain
- Office workers who are confined to a desk or computer for long periods of time
- People suffering from nervousness or insomnia

PACKAGE CONTENT


This package contains: the Main Unit, a User Manual and a Power Charger.


Main Unit

User Manual

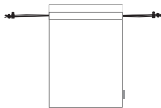Carrying Pouch

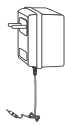Power Charger



Specifications are subject to change without prior notice. Refer to www.germanpool.com for the most up-to-date version.

PRODUCT STRUCTURE

Main Unit

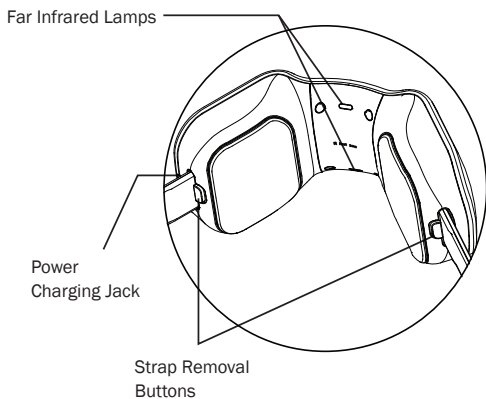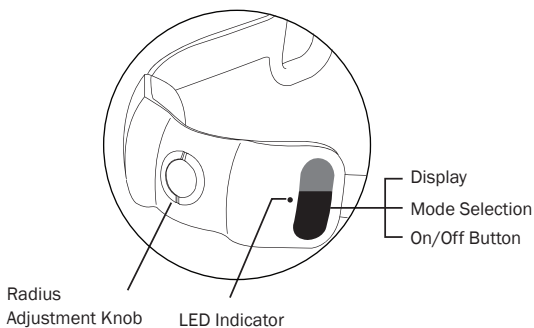Specifications are subject to change without prior notice. Refer to www.germanpool.com for the most up-to-date version.

OPERATING INSTRUCTIONS

Check before use

1. The battery has sufficient charge
2. The Power Charger used DOES indeed come with this product


Charging

1. Insert the Power Charger's DC cylindrical plug into the Power Charging Jack of the Main Unit. Connect the other end with an AC power source.
2. Charging is in process when the LED Indicator lights up. During charging, the battery icon on the Display will blink alternately. When the Neck Massager is fully charged, the blinking stops, and the battery icon appears as a full bar.
3. Charging takes approximately 5-6 hours.


Note: The product will not operate while it is connected to the Power Charger.
Remove the cable from the Power Charging Jack before starting the unit.

How to wear

Pull up the arrow marking on the Radius Adjustment Knob, and twist the protruding "V" to adjust the distance between the left and right massager heads.

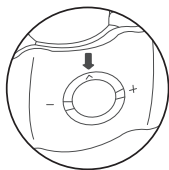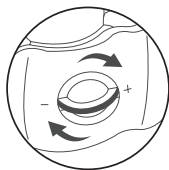Pluck and pull

Turn and adjust


Position the main unit directly on neck, hold strap handles firmly with both hands, and adjust until massager is in desired position. Adjust massage pressure and angle according to personal preferences.

Start

Press the On/Off  button. The product will start working, and the LED Indicator will display a red colour. The default massage cycle is 15 minutes, and the default massage mode is 1 (massage heads circle backward in tandem). To change massage mode, press Mode Selection button . Press the On/Off  button to select between massage cycles of 5, 10, and 15 minutes. The minutes left in the massage cycle will be shown on the Display.

The corresponding time and massage modes will be shown on the Display.

OPERATING INSTRUCTIONS

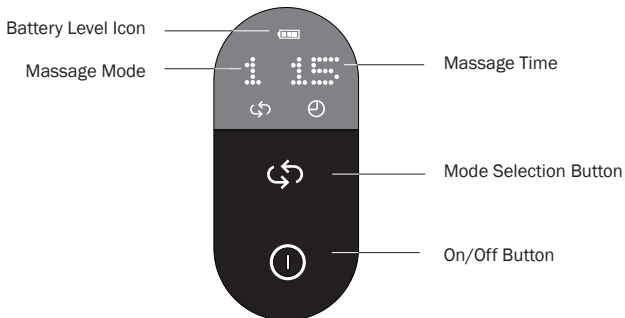When message cycle begins, the remaining time left in the message cycle will be shown on the Display. When message cycle is over the unit will automatically turn off.

Power-Off

To turn off the unit, long press the On/Off  button. The product will stop running and the LED Indicator will turn off.

Message Modes

The Neck Massager has four message modes. Select between them by pressing the mode selection  button.

Mode	Function
1	Message heads rotate backward in tandem
2	Message heads rotate forward in tandem
3	Both message heads rotate clockwise
4	Both message heads rotate anti-clockwise

CARE & MAINTENANCE

Storage

Place the unit in a safe place after use. Do not put heavy objects on top of the product. In case the storage period exceeds 3 months, please turn on the unit to check its function. Do not store the product close to a wall, on the floor, or on the roof. Keep in good ventilation conditions, away from direct sunshine and/or corrosive gases. Take measures to prevent damage to the unit from moisture, rats, and / or pollution.

Transportation

The product has a simple shockproof design and can be transported by air, railway, highway or ship. Avoid rain, snow, heavy impacts, and always keep the unit the right way up.

Check the power supply

Check the power supply before use. Do not switch on the device if the voltage exceeds the rated range of the device.

Maintenance

Contact a qualified technician if there is any fault with the unit. Do not disassemble or attempt to repair the device yourself.

TROUBLESHOOTING

If there is any problem when using the device, please follow the instructions below to check. If the problem persists or the device cannot function normally, please turn off the power supply and contact your local service centre. Do not disassemble or repair the device yourself, as risk of injury may occur.

Problem	Possible Cause	Solution
The product does not turn on	Battery is low	Charge at once
Massage strength is low	Battery is low	Charge at once

SPECIFICATIONS

Model	MN51
Dimensions	(H) 120mm (W)105mm (L) 220mm
Net Weight	855g
Power	5W
Voltage	DC 5V/1A

Specifications are subject to change without prior notice. Refer to www.germanpool.com for the most up-to-date version.

If there is any inconsistency or ambiguity between the English and Chinese versions, the Chinese version shall prevail.

請填寫以下空格並與發票一併保存，以便維修時核對資料。用戶提供本公司之所有資料必須正確無誤，否則免費保用將會無效。
Please fill out the form below. This information and the purchase invoice will be required for any repairs. Warranty will be invalid if information provided is found to be inaccurate.

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Warranty Terms & Conditions

Under normal operations, the product will receive a 1-year full warranty provided by German Pool (effective starting from the date of purchase).

* This warranty is not valid until Customer's registration information is returned to and received by our Service Centre within 10 days of purchase via one of the following methods:

- 1) Visit our website and register online: www.germanpool.com/warranty.
- 2) Complete all fields on the attached Warranty Card and mail it back to our Customer Service Centre along with a photocopy of the purchase invoice (Warranty Card and invoice copy will not be returned to Customer).

1. Customer who fails to present original purchase invoice will be disqualified for free warranty.
2. Customer should always follow the operating instructions. This warranty does not apply to:
 - labour costs for on-site installation, check-up, repair, replacement of parts and other transportation costs;
 - damages caused by accidents of any kind (including transport and others);
 - operating failures resulted from incorrect voltage, improper operations and unauthorized installations or repairs.
3. This warranty is invalid if:
 - purchase invoice is modified by unauthorized party;
 - product is used for any commercial or industrial applications;
 - product is repaired or modified by unauthorized personnel, or unauthorized parts are installed;
 - serial number is modified, damaged or removed from the product.
4. German Pool will, in its discretion, repair or replace any defective part.
5. This warranty will be void if there is any transfer of ownership from the original purchaser.

Customer Service & Repair Centre

German Pool (Hong Kong) Limited
8/F, China Travel (Cargo) Logistics Centre
1 Cheong Tung Road
Hung Hom, Kowloon, Hong Kong
Tel : +852 2333 6249
Fax : +852 2356 9798
Email : repairs@germanpool.com

此保用只適用於香港 | This warranty is valid only in Hong Kong

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Room 113, Newport Centre Phase II, 116 Ma Tau

Kok Road, Tokwawan, Kowloon, Hong Kong.

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© All rights reserved. Copying, reproducing or using the contents of this manual is not allowed without prior authorization from German Pool, violators will be prosecuted.

www.germanpool.com